AMT Bridle™

Nasal Tube Retaining System


An ISO 13485 Registered Company

Manufactured by:
Applied Medical Technology, Inc.
8000 Katherine Boulevard
Brecksville, OH 44141 USA

Phone: 1-440-717-4000
Toll Free: 1-800-869-7382
Fax: 1-440-717-4200

www.appliedmedical.net

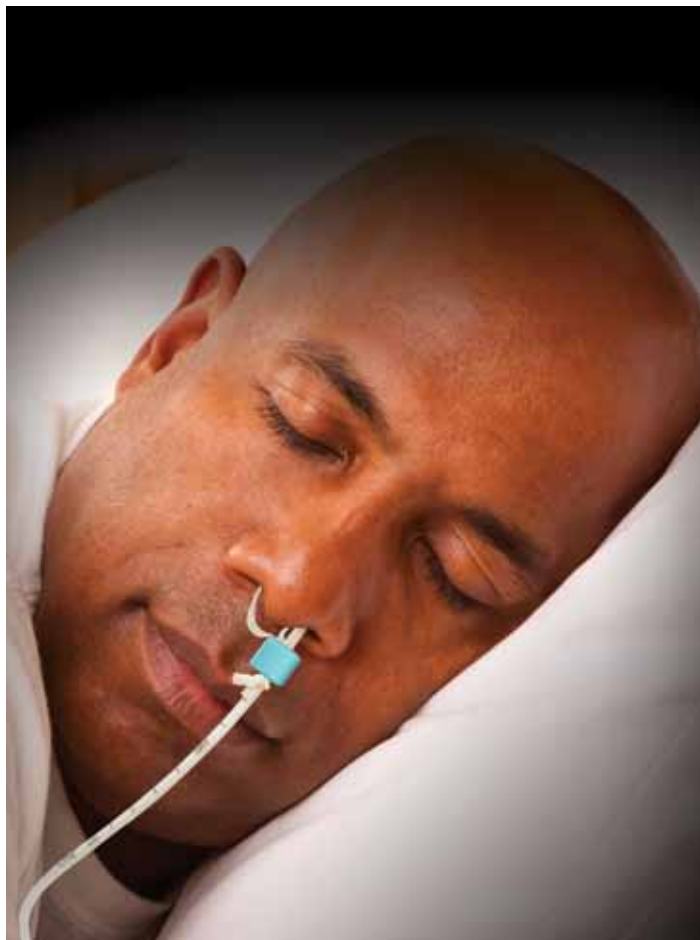AMT Bridle™ Nasal Tube Retaining System

ES	Sistema de Retención de Tubo Nasal, AMT Bridle™
DE	AMT Bridle™ Nasensonden-Haltesystem
IT	AMT Bridle™ Sistema di sostegno per tubo nasale
FR	Système de retenue pour tube nasal AMT Bridle™
NO	AMT Bridle™ System for tilbakeholding av Nesesonde
SV	AMT Bridle™ Fasthållningssystem för nästub

MR Safe


ES	Seguro para Resonancia Magnética
DE	MR-sicher
IT	Sicura per la RM
FR	Résistant aux champs magnétiques
NO	MR-sikkert
SV	MR-säker

DEHP Free Formulation

ES	Fórmula libre de DEHP (di(2-etilhexil) ftalato)
DE	DEHP (Di(2-ethylhexyl) Phthalat)-freie Formel
IT	Formulazione priva di DEHP (di(2-etilesil) ftalato)
FR	Formule sans DEHP (Phtalate de de(2-éthylhexyle))
NO	Ingen DEHP (Di(2-ethylhexyl) ftalaat)
SV	Sammansättning fri från DEHP (di(2-etylhexyl)ftalato)

Latex Free (No Latex)

ES	Libre de látex (sin látex)
DE	Latexfrei (kein Latex)
IT	Senza lattice
FR	Sans latex
NO	Lateksfri
SV	Latex-fri (Ingen Latex)

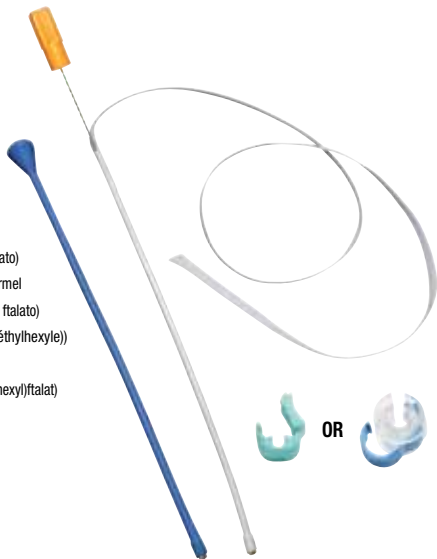AMT Bridle™ Nasal Tube Retaining System

Caution: Federal (USA) law restricts this device to sale, distribution, and use by or on the order of a physician.

A thorough nasal exam is recommended before AMT Bridle placement in patients to ensure that adequate bone remains in the posterior septum to support the AMT Bridle System.

A nasal exam should also be conducted after AMT Bridle placement to ensure that the umbilical tape loop has passed behind the nasal septum rather than through a septal perforation.

NOTE: System is supplied non-sterile, for single use only. Please inspect all contents of the kit for damage. If damaged, do not use product.

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INDICATIONS FOR USE

For use with Nasal Gastric and Nasal Enteral Feeding Tubes. To prevent inadvertent displacement or removal of nasal tube from patients requiring one.

CONTRAINDICATIONS FOR USE

This device is contraindicated for patients with nasal airway obstructions and abnormalities, and facial and/or cranial fractures. Do not use on patients that may pull on the AMT Bridle to such a degree as to cause serious injury.

Warning: Excessive traction on the AMT Bridle may cause tube displacement or nasal injury; an alternative means of securing the feeding tube should be explored.

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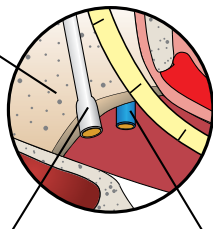
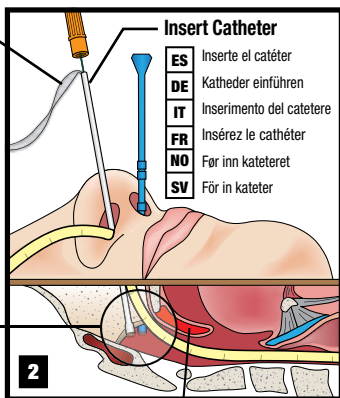
- 2** Insert the white Bridle catheter into the opposite nostril to approximate the magnets. Pull back on the orange stylet about 1 cm (1/2") until the magnets connect. The "click" of the magnets may be heard or felt. Equal lengths of both probes (minus the orange stylet) should be exposed.

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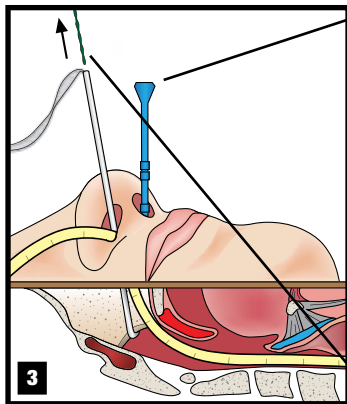
It may be necessary to advance or manipulate the probe to achieve magnet contact.

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Remove stylet AFTER magnet connection

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- 3** If necessary, gently twist the probes from side to side and/or up and down to encourage contact between the magnets. If no contact has occurred, then advance both the white catheter and the blue probe. **Important: Once contact has occurred, remove the orange stylet completely from the white catheter and dispose.**

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
Catheter and umbilical tape will advance to form loop

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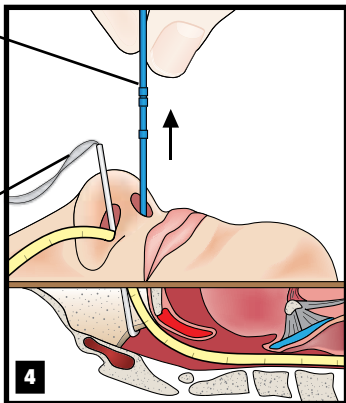- 4** Slowly withdraw the blue probe and allow the white catheter to advance through the nose. Continue until only the cloth umbilical tape is in the nose. This creates a loop or “bridle” around the vomer bone. If catheter and umbilical tape do not advance out of the opposite nostril, remove catheter, replace the stylet and start over at Step 1.

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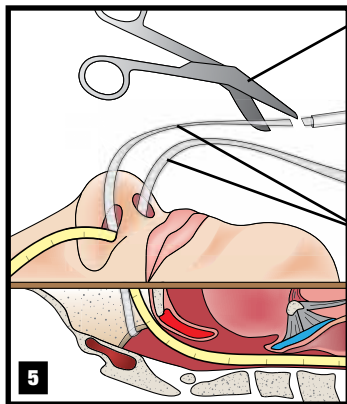
- 5** Cut the catheter portion from the umbilical tape. Dispose of white catheter and blue probe. If the nasal tube has not been placed, insert now.

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
- 6 Important:** Use caution when handling the AMT Bridle clip near patient's mouth; consider covering patient's mouth. Clip should be located close to nostril, approximately 1cm (1/2") or one "finger width" from nostril. **NOTE:** The clip should not touch the nostril. For clip sizes 5-12F: Place the nasal tube into the clip's channel or groove. **For clip sizes 14-18F refer to page 25.**

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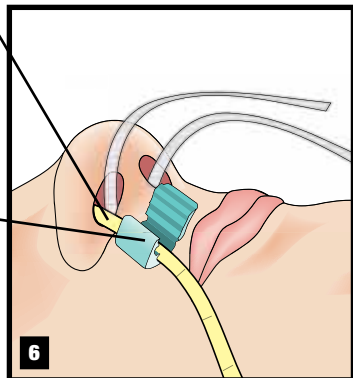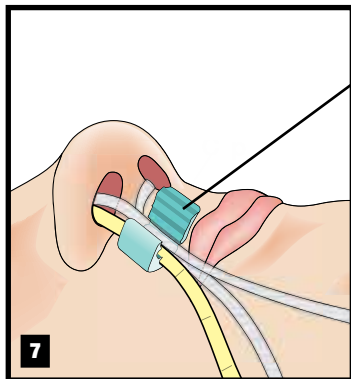
- 7** Ensure nasal tube fits securely in channel. Keeping them flat, place both strands of umbilical tape in the clip's hinge.

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
- 8** Close the clip by folding over the plastic edge and firmly snapping it shut.

NOTE: The clip is difficult to re-open, so ensure proper position of the nasal tube, umbilical tape and clip prior to closure.

ALSO NOTE: Some nasal tubes are more rigid than others and may require a higher force to close clip. **WARNING:** Do not use any foreign object to open or close clip as this may damage the clip, making it less secure and/or ineffective. If damaged, use a new clip.

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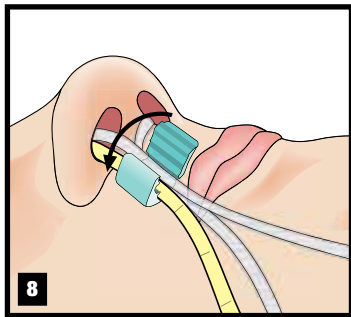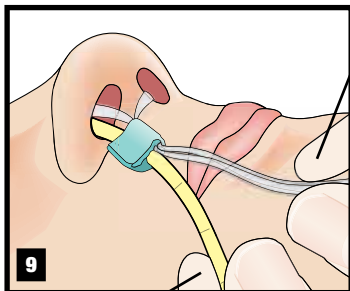
Gently pull on umbilical tape ends to verify clip is closed.

- ES** Hale suavemente los extremos de la cinta umbilical para verificar que el clip está cerrado.
- DE** Vorsichtig an den Enden der Sicherungsbänder ziehen, um die Funktionsfähigkeit der geschlossenen Klemme sicherzustellen.
- IT** Tirare con cautela le estremità del nastro ombelicale per verificare che la clip sia chiusa.
- FR** Tirez doucement sur les extrémités du pansement ombilical pour vérifier que l'agrafe est bien fermée.
- NO** Dra forsiktig på endene til den medisinske tapen for å kontrollere at klemmen er lukket.
- SV** Dra sakt på navelsträngstejpen ändrar för att bekräfta att klämman är stängd.

Hold tube fixed.


- ES** Mantenga fijo el tubo.
- DE** Die Nasensonde festhalten.
- IT** Mantenere fermo il tubo.
- FR** Retenez le tube fixé.
- NO** Hold sonden fast.
- SV** Håll tuben fixerad.

9 After the clip has been fully closed, tie the two umbilical tape ends together (excluding the tube) creating a simple knot. Repeat 2-3 times.

- ES** 9. Después que el clip está completamente cerrado, amarre los dos extremos de la cinta umbilical para unirlos (excluyendo el tubo) creando un nudo sencillo. Repita 2 a 3 veces.
- DE** 9. Sobald die Klemme fest geschlossen ist, die beiden Sicherungsbänderenden miteinander verknöten (nicht mit Nasensonde). Insgesamt 2-3 Knoten machen.
- IT** 9. Legare le due estremità del nastro ombelicale insieme (escludendo il tubo), dopo che la clip è stata chiusa completamente, effettuando un semplice nodo. Ripetere per 2 - 3 volte.
- FR** 9. Une fois que l'agrafe a été complètement fermée, attachez les deux extrémités du pansement ombilical ensemble (sans prendre le tube) en faisant un nœud simple. Nouez ainsi 2-3 fois.
- NO** 9. Etter at klemmen har blitt lukket skikkelig, knyt sammen de to endene i den medisinske tapen (bortsett fra sonden) ved å lage en enkel knute. Gjenta 2-3 ganger.
- SV** 9. Efter att klämman helt har stängts, knyt ihop de två ändarna på navelsträngstejpen (exkl. tuben) och skapa en enkel knut. Upprepa 2-3 gånger.

Cut excess umbilical tape.

- ES** Corte la cinta umbilical sobrante.
- DE** Überschüssiges Sicherungsband in der Länge nach Bedarf abschneiden.
- IT** Tagliare il nastro ombelicale in eccesso.
- FR** Coupez l'excédent de pansement ombilical.
- NO** Klipp av overflødig medisinsk tape.
- SV** Skär av överflödig navelsträngstejp.

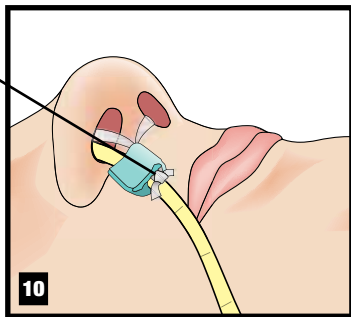- 10** The AMT Bridle™ may remain in place for the entire time the feeding tube is in place.

Note the placement of the AMT Bridle in the patient's chart. **WARNINGS:** Throughout the useful life of the AMT Bridle, the clip must be visually inspected for signs of damage or loose attachment to the nasal tube. Monitor the clip and nasal tube position regularly for feeding tube migration.

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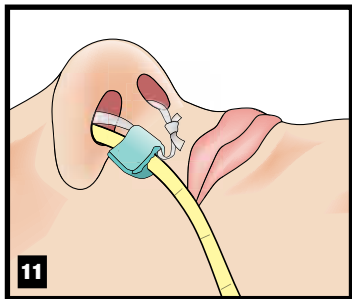
Optional Method

11 As an option, a single strand of umbilical tape may be placed in the clip and tied off.

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12 Clip Opening Tool

If the clip must be opened, place the side of the opening tool into the clip. Push in and turn slightly to open.

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
13 Removal

To remove the AMT Bridle™ and tube:
Cut only ONE strand of the umbilical tape.
Gently pull both the Bridle and nasal tube
out of the nose. To remove only AMT Bridle:
Cut both ends of umbilical tape and remove
(clip remains in place). **CAUTION: Great
care must be taken to prevent accidental
swallowing as the tape can easily slide
back into the nose.**

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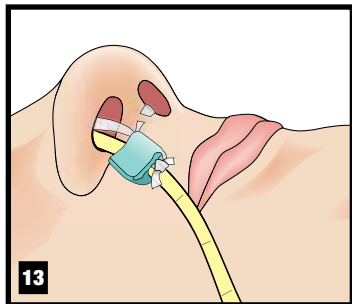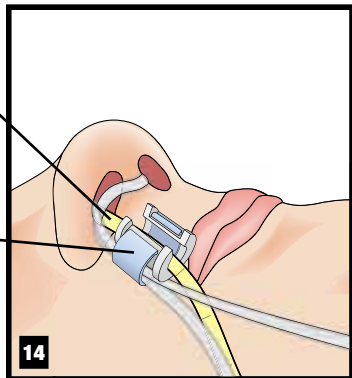
For Clip Sizes 14–18F

Ensure nasal tube fits securely in channel. Keeping them flat, place both strands of umbilical tape between the soft inner portion of the clip and the outer, more rigid plastic section. Go back to Step 8.

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Manufactured by:
Applied Medical Technology, Inc.
8000 Katherine Boulevard
Brecksville, OH 44141 USA
Phone: 1-440-717-4000
Toll Free: 1-800-869-7382
Fax: 1-440-717-4200

www.appliedmedical.net

C-4281-H 12/12